


# Alpha® Swab with Long Handle

## Alpha Series

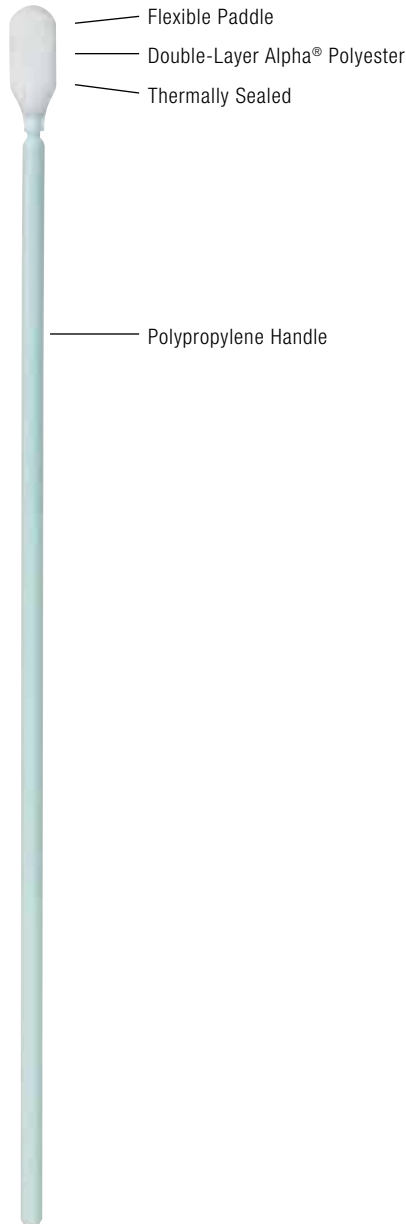

Actual Size

### Description

The Alpha® Swab with Long Handle is an excellent general purpose swab for cleaning surfaces and hard-to-reach areas. Its laundered, knitted Alpha® polyester head is extremely clean. A long handle and flexible internal head paddle provide good support. This swab can replace typical long-handled cotton swabs in applications requiring superior cleanliness. The Alpha® Swab with Long Handle is also widely recommended for surface sampling and cleaning-validation work. Its cleanliness provides a neutral background for the sample obtained.

### Features

- Cleanroom-laundered, double-layer knitted Alpha® polyester tip over flexible paddle
- Long, light green polypropylene handle
- Thermally bonded head
- Soft and nonabrasive
- Excellent chemical resistance

### Benefits

- Ultralow nonvolatile residue
- Low particle and fiber generation
- Good sorbency
- No contaminating adhesives
- Compatible with many solvents
- May be sterilized by autoclave

### Applications

- Cleaning with solvents such as IPA
- Sampling surfaces for TOC or cleaning validation
- Applying lubricants and other liquids
- Removing excess materials

### Products

TX Number	Description	Packaging
<b>TX761</b>	Alpha® Swab with Long Handle	Bag-Within-A-Bag® 100 swabs/bag 1 inner bag of 100 swabs 10 bags/case



**North America**  
1210 South Park Drive  
Kemersville, NC 27284  
Tel 336.996.7046  
Fax 336.996.2297  
[www.texwipe.com](http://www.texwipe.com)  
[info@texwipe.com](mailto:info@texwipe.com)

**Europe/Middle East**  
Saffierlaan 5  
2132 VZ Hoofddorp  
The Netherlands  
Tel +31 88 1307 400  
Fax +31 88 1307 499

**Asia/Pacific**  
50 Tagore Lane  
#02-01 Entrepreneur Centre  
Singapore 787494  
Tel +65 6468 9433  
Fax +65 6468 6772

# Alpha<sup>®</sup> Swab with Long Handle

## Physical Characteristics

Head material	Knitted Alpha <sup>®</sup> polyester
Head width	6.8 mm (0.268")
Head thickness	2.8 mm (0.110")
Head length	16.8 mm (0.661")
Handle material	Polypropylene
Handle width	3.2 mm (0.126")
Handle thickness	3.2 mm (0.126")
Handle length	145.5 mm (5.728")
Total swab length	162.3 mm (6.390")
Head bond	Thermal
Handle color	Light green
Design notes	Flexible head paddle; long handle

## Contamination Characteristics

Property	Measurement	Test Method*
Nonvolatile residue		
IPA extractant	0.02 mg/swab	
Ions		
Chloride	0.01 µg/swab	
Sulfate	< 0.01 µg/swab	
Nitrate	< 0.01 µg/swab	
Phosphate	< 0.01 µg/swab	TM2: Laboratory Testing For Swabs
Fluoride	< 0.01 µg/swab	
Potassium	0.01 µg/swab	
Calcium	< 0.01 µg/swab	
Sodium	0.01 µg/swab	
Magnesium	0.01 µg/swab	

\*Test procedures available upon request.